

umbrella

Mobile Personal Emergency
Response System
mPERS

Getting Started






Initial Setup

1

Carefully take the device, the charging cradle, its cable, and the power supply out of the box.

2

Connect the charging cable to the cradle and to the power supply. Then connect the power supply to the outlet.

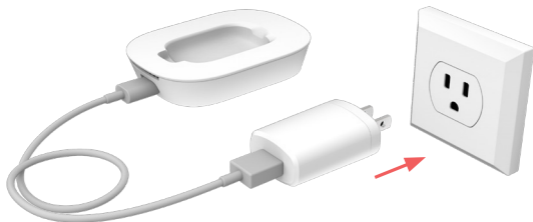

3

Place the device in the cradle for at least six hours.



The device will power on while charging. This may take some time.

4

Perform an initial test:

Pick up the device. Press and hold the panic button until you hear the message “emergency call”.

Speak to the monitoring service and confirm your location.



To initiate a call, **press and hold the panic button** until you hear the message “emergency call”.



To cancel or hang up, press twice

In an emergency, do not cancel or hang up the call before speaking with the monitoring service.

Caution- terminating an emergency call will result in the monitoring service not being informed of an emergency.

To check battery and reception levels.





Hold for 3 seconds

Check Umbrella device battery level and reception regularly. Charge nightly or as soon as possible if the battery level is low. A depleted battery or low reception may prevent your device from operating and communicating with the monitoring service, including communicating your emergency. When the battery level is low place the device in the cradle, as shown in step 3 of the initial setup.



Day-to-Day Use



Wear around your neck or on your belt.

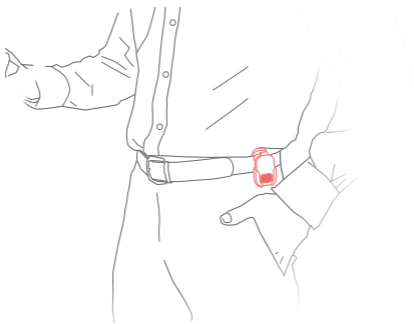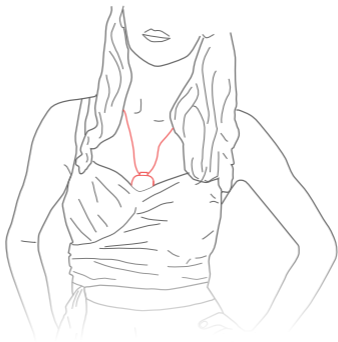



Wearing Umbrella on a Belt

1 Remove the cord.

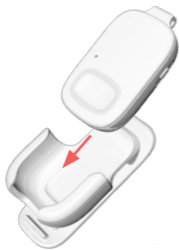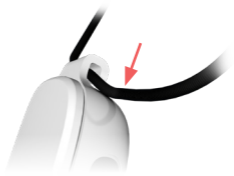

2 Slide Umbrella into the belt clip.



3 Clip onto belt or pocket.

For more information, see www.essencesmartcare.com/assistance-center

www.essencesmartcare.com • info@essencesmartcare.com



© 2024 All rights reserved to Essence SmartCare Ltd.



ESUC05259 Rev. 1.0

